





**1** Place the slide inside pastel

**2** Close the door

**3** After one minute your slide is ready to interpret

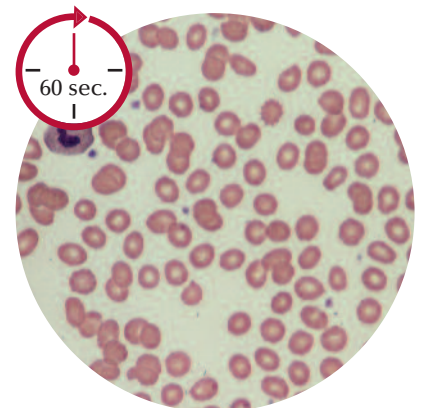

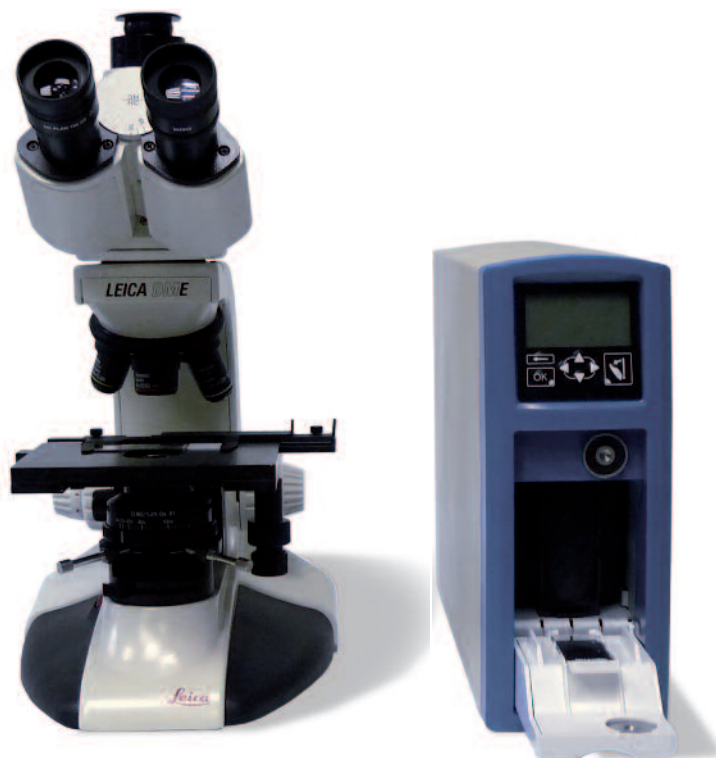

## VERY COMPACT INSTRUMENT

### READY TO USE

- Cartridge of reagents inside the stainer

### NO WASTE OF REAGENTS

- Storage of the reagents in darkness condition
- No contamination between reagents during the cycle
- Single use of the reagents

### CLEAN AND ENVIRONMENTAL SAFETY

- Waste integrated in the cartridge
- No toxic evaporation

### FAST

- 1 slide in 1 minute (include drying)

### UP TO 9 DIFFERENT STAINING CYCLES, AMONG THEM MALARIA CYCLE

### USER'S MAINTENANCE FREE

## PASTEL SPECIFICATIONS

Size (w x h x l) in mm	105 x 270 x 380
Weight (in Kg)	3
Power	External power supply 12 volts
Time of cycle	60 sec.
Consumable	Cartridge
Display	Graphic display

## CARTRIDGE SPECIFICATIONS

Size (w x h x l) in mm	70 x 80 x 220
Number of tests	50
Dyeing solution	RAL 555
Identification	RFID (Radio Frequency Identification)
Expiry date	2 years